
**NAIPO**  
FIND WELL-BEING

## User Manual

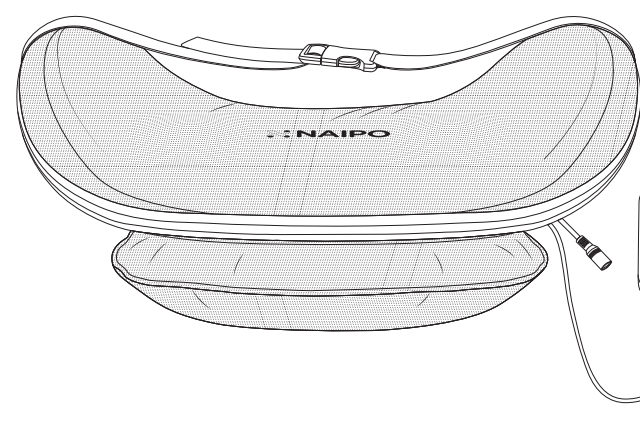**Lumbar Traction Massager with Vibration and Heat**  
Model: MGBK-Q1

Shenzhen Pinlong E-Commerce Co., Ltd.  
Web: www.naipocare.com  
E-mail: support@naipocare.com  
Add: Huanan City Electronic Trading Building P07  
Room 120, Longgang District, Shenzhen, China.

## Contents

English	01-03
Deutsch	04-06
Français	07-09
Español	10-12
Italiano	13-15

### Introduction

Thank you for purchasing the Nairo Lumbar Traction Massager with Vibration and Heat. Please read this user manual carefully and keep it in a safe place for future reference. If you need any assistance, please contact our support team with your product model number and Amazon order number.

### Included

- 1 x Lumbar traction massager
- 1 x User manual
- 1 x Adapter

### Product Specifications

Model	MGBK-Q1
Name	Nairo Lumbar Traction Massager with Vibration and Heat
Input voltage	12V 1A
Power consumption	12W
Timer	15 minutes
Load capacity	100KG
Surface heating temperature	≤45°C


### Important Safeguards

- People with back diseases should consult a doctor before using the product. When using the product for the first time, we recommend you manually inflate or deflate the air bag, and 50% of full inflation is enough. When using the product for the first time, please do not inflate the air bag excessively to avoid hurting the waist.
- Do not use the product for more than 30 minutes at a time.
- Heating function: It takes about 1 minute for the product to warm up.

### Q & A


Problem	Solution
The product does not work / The product stops working.	<p>Please check whether the adapter is connected tightly to the power supply.</p> <p>Please check if the adapter is connected to the product properly.</p> <p>The product shuts off after 15 minutes of consistent use. Please restart the product.</p> <p>Overheat protection is activated if the product has been in use for long (especially when the product has been used continuously for several sessions, with 15 minutes for each session). Please disconnect the product from the power supply and let it cool down for about 1 hour.</p> <p>Please contact the customer service if the above methods fail to solve the problems.</p>

### Control Panel

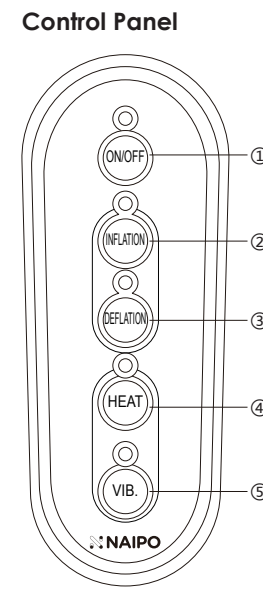
- 1. Power Button**  
After connecting the product to the power supply, the indicator near the power button will flicker (red), indicating that the product is under the standby mode. Press the power button to turn on the product with the indicator on all along, and the product will shift to automatic working mode, under which inflation, deflation, heat, and vibration would work automatically.
- 2. Inflation Button**  
Press the button to shift to manual mode. Long press the button to inflate the air bag consistently. After the air bag has been fully inflated, the inflation will stop.
- 3. Deflation Button**  
Press the button to shift to manual mode. Long press the button to deflate the air bag consistently. After the air bag has been completely deflated, the deflation will stop.
- 4. Heating Button**  
When the product is turned on, the heating function will be activated by default with the indicator near the button on (red); Press the button to turn off the heating function with the indicator off.
- 5. Vibration Button**  
When the product is turned on, the vibrating function will be activated by default with the indicator near the button on (red); Press the button to turn off the vibrating function with the indicator off.


### Product Diagram
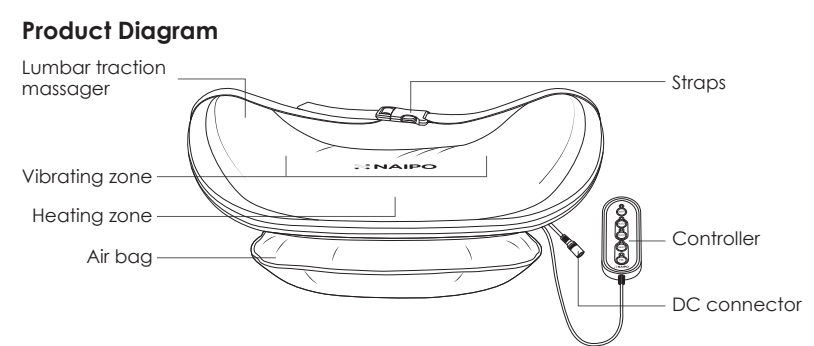
### Product Operation Guide

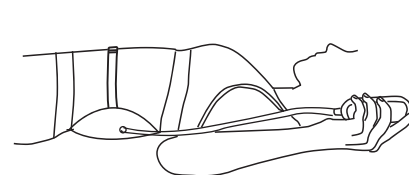- Fasten and stabilize with the straps, and then lie down slowly. After lying and adjusting the position, begin the massage and control the working modes via the controller.

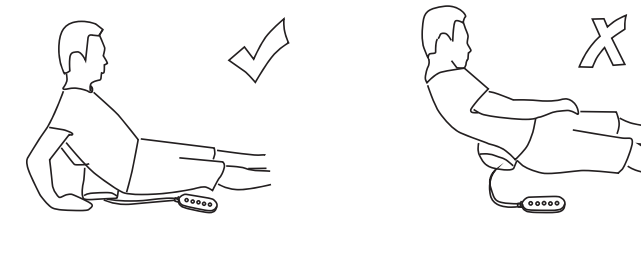- The product will shut off automatically after 15 minutes of massage, or you can stop the massage in advance by pressing the power button. After the massage has finished and the air bag has been completely deflated, unfasten the straps and continue to lie for 1-3 minutes before getting up slowly with support of arms.

### Product Maintenance

- Water or other liquids should never be let in the product, otherwise malfunction or damage might occur.
- Wipe the product with soft and mildly damp cloth for cleaning.
- Please store the product in a clean, dry, and cool place. Avoid exposing the product to sunlight for long.

### Warranty & Customer Support

For questions, assistance, support, or warranty claims, please contact us via the email or the address below.  
Web: www.naipocare.com  
E-mail: support@naipocare.com  
Address: Huanan City Electronic Trading Building P07 Room 120, Longgang District, Shenzhen, China.

\*Note: Products purchased directly from Nairo benefit from our after-sales service. If the purchases are made from other selling teams, we recommend you contact them for warranty issues.

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