
Shenzhen Bestvap Technology Co., Ltd.  
**Web:** www.naipocore.com  
**E-mail:** bestvap@naipocore.com  
**Add:** 503 Mantao Building 720 Minzhi Avenue Minzhi Street Longhua District, Shenzhen, Guangdong 518006

**CE** **RoHS**

**NAIPO**  
 FIND WELL-BEING

**User Manual**

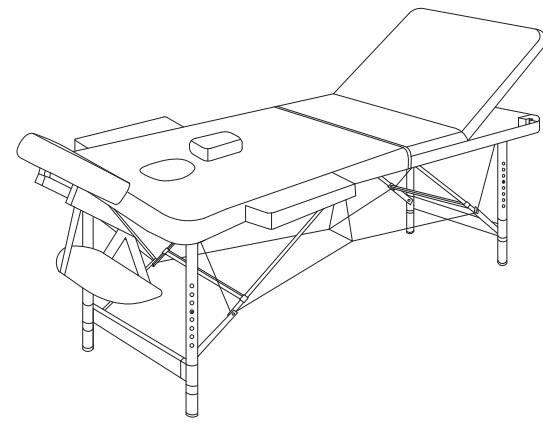**Three-fold Portable Massage Table with Aluminum Feet**

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
**Introduction**  
 Thank you for purchasing the Naipo Three-fold Portable Massage Table with Aluminum Feet. Please read this user manual carefully and keep it in a safe place for future reference. If you need any assistance, please contact our support team with your product model number and Amazon order number.

**Included:**  
 1 x Massage table  
 1 x Face hole plug  
 1 x Face cradle cushion  
 1 x Face cradle  
 1 x Armrest  
 2 x Side armrest  
 1 x Carrying bag  
 1 x User manual

**Product Specifications**

Model	MG8C-301H
Name	Naipo Three-fold Portable Massage Table with Aluminum Feet
Dimensions	185*70*(60-81)cm (adjustable height)
Maximum working weight capacity	270kg


**Using Guide**

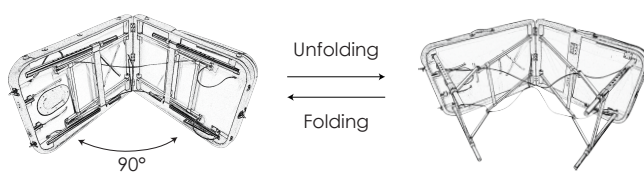- 01 -

**User tips:**  
 a. When unfolding, lay the folded product on its side, unfold it until the angle is 90°, and then unfold the aluminum feet.  
 b. When folding, lay the unfolded product on its side, fold it until the angle is 90°, and then fold the aluminum feet.  
 c. Please make sure that the wire ropes do not wind around the aluminum feet when unfolding or folding the product.


**Installation of Face Cradle Cushion**

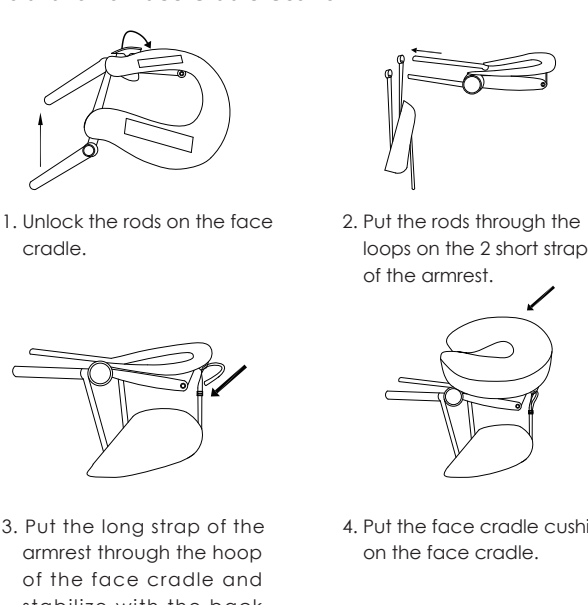
1. Unlock the rods on the face cradle.  
 2. Put the rods through the loops on the 2 short straps of the armrest.  
 3. Put the long strap of the armrest through the hoop of the face cradle and stabilize with the back buckle.  
 4. Put the face cradle cushion on the face cradle.

- 02 -

**Height Adjustment**  
**Instructions:**  
 After pressing the button on the aluminum foot into the small holes, grip the upper part of the aluminum foot with one hand, and move the lower part of the aluminum foot with another hand to desired height. When the height is selected, let the button stay on the corresponding hole for stability.

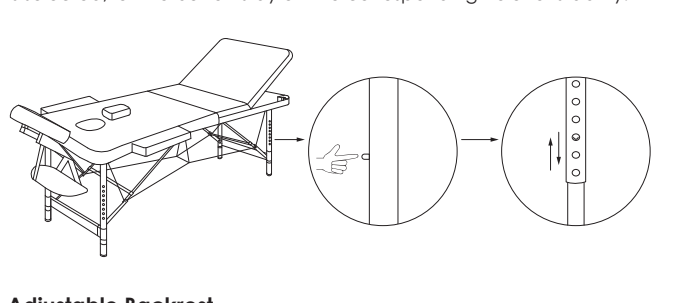**Adjustable Backrest**  
**Lifting:** Gently lift this part with 2 hands, and the lifter would respond accordingly, which allows the selection of desired height.  
**Lowering:** Lift this part to the maximum height and then put it down.

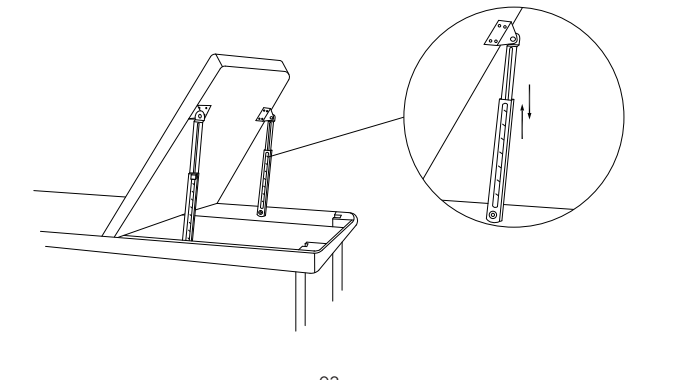- 03 -

**Product Maintenance**  
 a. The product is designed for massage, not for intense activities.  
 b. Please wipe the surface of the product with wet sponge or mildly damp cloth for cleaning.  
 c. Please do not clean the surface of the product with any corrosive detergent.  
 d. Please fold and store the product if it's not to be used for long.

**Warranty & Customer Support**  
 If you need help regarding product issues, after-sales service, or other relevant problems, please contact us via the email on the cover.

\* Note: Products purchased directly from Naipo will enjoy our after-sales service. If the purchases are made from other selling teams, we recommend you contact them for warranty issues.

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