


正面

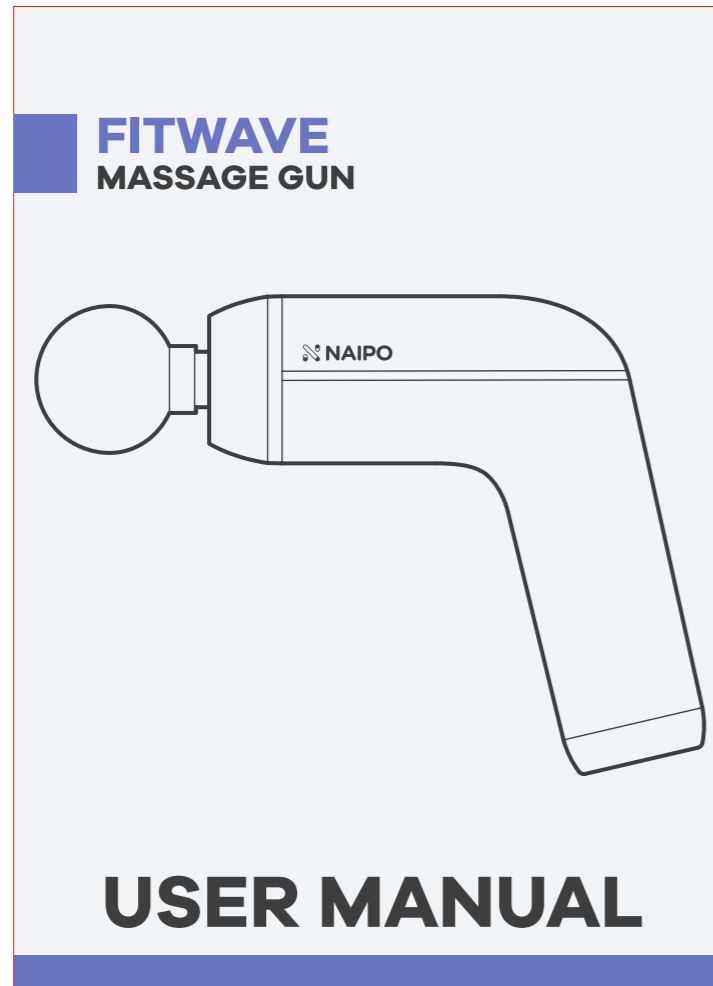

CONTENT

FITWAVE MASSAGE GUN

<p>02 INTRODUCTION</p> <p>04 MESSAGE AREA</p> <p>05 OPERATION GUIDE</p> <p>06 MALFUNCTIONING AND MAINTENANCE</p> <p>08 LEGAL NOTES</p>	<p>02 DANGER & WARNING</p> <p>04 IN THE BOX</p> <p>06 CLEANING AND STORAGE</p> <p>07 ENVIRONMENTAL PROTECTION</p> <p>08 SPECIFICATIONS</p>
---	---

INTRODUCTION

Thank you for purchasing NAIPO FitWave Massage Gun. This product is suitable for a physical massage before or after exercise. It can reduce muscle stiffness, relieve muscle soreness, increase blood circulation, relieve fatigue and promote muscle recovery.

Before use, please read this manual carefully to ensure correct operation. Please keep this manual for future reference.

DANGER & WARNING

Attention - Read all instructions carefully before using the product. Failure to follow the warnings and instructions may cause electric shock, fire or personal injury.

A) Safe massage environment

- Do not operate the NAIPO FitWave Massage Gun in a potentially explosive environment, such as around flammable liquid, gas, or dust.
- During operation, please keep away from children; any distraction may cause you to lose control.

B) Personal safety

- Do not use the massager when you are tired or under the influence of drugs or alcohol. If you are negligent when operating the massage equipment, it may cause personal injury.
- Always maintain proper body balance, so as to better control the massage equipment in unexpected situations.

ATTENTION - READ ALL INSTRUCTIONS CAREFULLY BEFORE USING THE PRODUCT:

- Do not use this product on the head or other bony body parts. Please use only on the soft tissue parts of the body. It is not recommended for use on areas with less muscle coverage.
- The product cannot replace medical treatment. If you have physical pain, are receiving medical treatment, or are using other medical devices, please consult your doctor before using this product.

- To avoid pinching, keep fingers and hair away from the back of operating massage head.
- Do not put the product in the water. If the product drops into water, please send it away for inspection.
- Keep this product away from fire.
- Do not disassemble this product in any way.
- Please switch the power off when the product is not in use.
- Do not use the product while driving, and observe all driving regulations to ensure safety.
- This appliance can be used by people with reduced physical, sensory or mental capabilities or lack of experience and knowledge if they have been given supervision or instruction concerning use of the appliance in a safe way and understand the hazards involved. Children should not play with the appliance. Cleaning and user maintenance should not be done by children without supervision.

When the massage head is not installed, do not touch this part so as to avoid pinching.
When installing the massage head to start the device, do not touch this part so as to avoid pinching.

PEOPLE NOT RECOMMENDED TO USE THE PRODUCT

If you have any of the following conditions, please consult a doctor before using this product:

- People who are pregnant;
- People with pacemakers or concerns about their health;
- People with physical disabilities;
- People who recently had surgeries.

These conditions do not mean that you cannot use the product for exercise, but we recommend that you consult a doctor first.

Practical experience has shown that in many cases, adding vibration training to a treatment plan can also be an effective treatment method. This must be done under the advice or supervision of a doctor, expert, or physical therapist.

RECOMMEND MESSAGE AREAS

DO NOT MESSAGE THESE AREAS

● Avoid massaging areas circled in red.

IN THE BOX

反面

CE FC RoHS

CONTACT US

Shenzhen Bestvap Technology Co., Ltd.
 Web: www.naipocare.com
 Email: bestvap@naipocare.com
 Add: 503 ManTao Building 720 Minzhi Avenue Minzhi Street Longhua District, Shenzhen, Guangdong 518000

MADE IN CHINA

OPERATION GUIDE

1 CHARGING
Connect to power supply with USB-C Cable until the battery is fully charged.
* Full charge is recommended for first use

2 MESSAGE HEAD INSTALLATION
Attach the desired massage head.
Round Head | Arms, Waist, Buttocks, Thighs, Calves
Y-Shaped Head | Neck, Either Side of Spine, Achilles Tendon
Bullet Head | Meridians, Joints, Palms, Soles
Pneumatic Head | Shoulder, Neck, Abdomen

3 MESSAGE
Turn On
Press and hold the start button for 2 seconds until the light is on.
Start Massage
Press again the start button to start.
Press the button to adjust the speed from Level 1 - 5.
Turn Off
Press and hold the start button for 2 seconds.

4 INTELLIGENT AUTO-ADAPT SPEED
Designed with Intelligent Auto-adapt Technology. FitWave Massage Gun automatically adjusts the intensity based on your usage.

CHARGING
1. Charge when the indicator red light is flashing 3 times, as it indicates the battery is nearly empty.
2. When charging, the indicator light flashes in white. When the battery is full, it shows 5 white lights.

ATTENTION: DO NOT USE THE DEVICE WHEN CHARGING.

CLEANING AND STORAGE

- Do not exceed the recommended daily usage of 30 minutes per day.
- Make sure the product has been turned off after use or before cleaning.
- Do not throw during use.
- Clean with a dry cloth, do not use engine oil, organic solvents, or other chemical agents to clean this product.
- Store the product in a dry place when not in use for a long time. Keep away from environments with high humidity or temperature. It is recommended to charge the product every six months to prevent damage to the battery due to power loss.

MALFUNCTIONING AND MAINTENANCE

PROBLEMS	POSSIBLE REASONS	SOLUTIONS
The product does not turn on	The temperature is too low: The product's battery activity will decrease and fail to start normally if it is in an environment below 0 degrees Celsius for a long time. The machine may appear able to start normally but cannot operate normally. The temperature is too high: The product may be overheated if it is used in an environment above 40 degrees Celsius. The machine may not be able to start normally. Insufficient battery power: When the battery is in a low power state and the machine is not charged for a long period of time, the battery may enter the low power protection (over-discharge protection). At this time, the machine will not work normally.	Put the device in an environment above 0 degrees Celsius for 2 hours, or use a hairdryer to quickly increase the machine's temperature. Put the device in an environment under 40 degrees Celsius for 2 hours, then turn it on. In the shutdown state, plug the device in to charge. The machine can work normally after being fully charged.
Massage heads fall off	Massage heads are not attached properly.	Check and attach the massage heads correctly.
The product will not charge & the panel light is off	The charging cable, charging port, and power supply are not correctly in place.	Check and reconnect the charging cable, charging port, and power supply properly.

Slight abnormalities in use	New machine breaking in period.	Use it for some time to observe whether there is still any abnormality. If so, please send it away for inspection.
Stops working after 15 minutes of use.	The product's 15-minute automatic shutdown timer function.	Restart the machine by pressing the start button.
Stops working when the massage head is blocked.	The product starts abnormal operation protection pattern and switches from operation to standby mode.	Restart the machine by pressing the start button.

ATTENTION: If the above methods do not solve the corresponding problem of the product, please contact NAIPO customer service.

Do not disassemble or modify the product by yourself. NAIPO will not be responsible for any loss caused by unauthorized disassembly or modification of the product. This action may put the user in danger and / or render the product unusable.

DISCLAIMER

NAIPO shall not be liable for costs, damages or repairs incurred as a result of, due to, or from:

- The warranty period having expired.
- Damage caused by force majeure, improper use or transportation.
- Unauthorized actions such as: product repair, misuse, abuse, modification, disassembly, etc.
- Do not follow the instructions.
- Damaged of the inner and outer packaging of the product, missing parts, etc.
- Inability to provide proof of purchase.
- External and internal wear caused by mechanical operation.

The warranty period is 12 months, and the authority is not transferable. The user can enjoy the warranty service with valid documents or certificates within the warranty period from the date of purchase.

ENVIRONMENTAL PROTECTION

This product contains batteries and must not be disposed of with household garbage. Please obtain recycling recommendations according to local environmental regulations.

LEGAL NOTES

All the above content does not contain any treatment statement. The diagnosis and treatment of the disease requires the help of a doctor or therapist. NAIPO does not promise that the product has any curative effect.

SPECIFICATIONS

PRODUCT NAME	NAIPO FitWave Massage Gun
MODEL	NPMG-2008
NET WEIGHT	1.34lb / 610g
PRODUCT DIMENSIONS	6.4×6.2×2.3in / 163×157×58mm
BATTERY CAPACITY	2600mAh
INPUT	5V = 2A
POWER	10W

虚线为折线，实线为切割线

成品尺寸：95x130mm
 材质：157g哑粉纸