FIND WELL-BEING

USER MANUAL

CONTACT US

Aukey International Limited

Web: www.naipocare.com

E-mail: wholesale@naipocare.com

Add.: Room 601 6/F Tung Hip Commercial Building,
248 Des Voeux Road Central Sheung Wan, Hong Kong

HQ: No. 102, Building P09, China South City, No. 1 Huanan Avenue,
Pinghu Street, Longgang District, Shenzhen, China

MADE IN CHINA

[/naipocare](#) [/naipocare](#) [/naipocare](#)

oPillow Plus

BACK MASSAGER

CONTENTS |

English

02

Deutsch

10

Français

18

Español

26

Italiano

34

IMPORTANT SAFETY INSTRUCTIONS

When using the product, all the following instructions and basic precautions should be taken.

DANGER

To reduce the risk of electric shock, burns, fire, etc.:
Always unplug this product from the electrical outlet immediately after use and before cleaning.
Please do not disassemble the product by yourself.

WARNING

Danger - To reduce risk of electric shock:

Do not operate under a blanket or a pillow; allow for normal heat dissipation so the product does not overheat, which could lead to fire, electric shock, or injury.

Do not turn off by pulling the power cord.

Only use this product with the power supply unit provided, and only with the voltage specified on the label.

Do not operate in places where corrosive products are used or where oxygen is administered.

Do not use in environments with high humidity (eg: bathroom) or high temperature, which could lead to electric shock or other safety hazards.

To avoid high temperatures, do not leave product in the car.

Do not use outdoors.

Attention - Read all instructions carefully before using the product:

This product is a personal massage product, not a medical instrument, and it cannot be used for other commercial activities.

Only use this product for the intended purpose as described in this manual.

Please switch the power off when the product is not in use.

The product should never be left unattended when plugged in. Unplug from outlet when not in use.

The product is equipped with a heating function. People insensitive to heat should be careful when using.

Keep the cord away from heated surfaces, so as not to cause serious burns.

Please pay special attention to high wear components such as the massage heads. Replace defective components immediately and/or keep the product out of use until repair.

Do not use the product if it has a damaged cord or plug, if it is not working normally, if it has been dropped or damaged, or if it has dropped into water.

Stop using the product immediately if you feel discomfort or pain during use.

Stop using the product if it is not working normally.

Do not use the product while driving, and observe all road regulations.

This product is not suitable for children, and children should not play with the product.

If you have any of the following conditions, please consult doctor before using this product:

People with physical disabilities.

Those who are developmentally disabled or have psychiatric issues.

Those with chronic illnesses.

Those who are pregnant.

Those who are menstruating.


Those with physiological issues.


People with pacemakers or concerns about their health.

MASSAGE AREAS

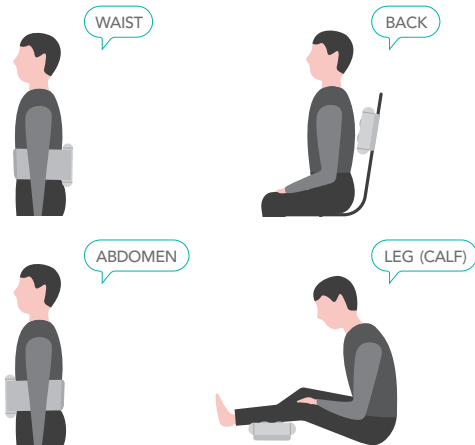IN THE BOX

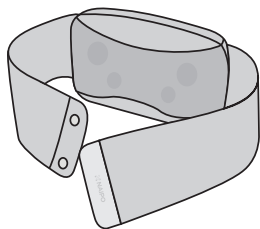Massager


Wall Adapter

OPERATION GUIDE

How to use

- Position the product on desired place.
- Hold down the power button for 2 seconds to start the product. The heating and rollers will automatically switch on. The default timer is 15 minutes.
- Press the heating button to turn off the heating function; (when the heating function is on, the indicator light is green).
- Press the speed adjustment button to choose the desired massage speed level: Low / Medium / High.
- Press the power button to stop the massage at any time.

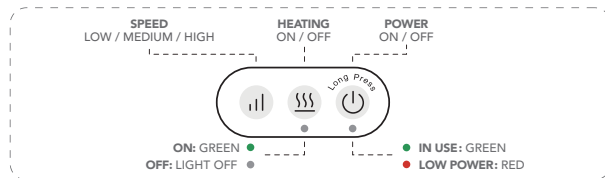Charging

- The indicator will light up in red when the product is low on battery. When the battery is low, plug in and charge.
- Insert the cylindrical end of the charger cable into the product charging port, and plug the adaptor into an electrical power outlet.
- A full charge cycle takes approximately 4 hours. During charging, the indicator flashes red and it turns to green when the product is fully charged.
- When connected to a power supply, the product will automatically charge if the battery is low.
- The product can be used continuously for 60 minutes in a fully-charged condition.

NOTE: The product can be used while plugged into a power source but the product will not charge.

CLEANING AND STORAGE

To prevent damage and discoloration, do not clean the product with corrosive detergent. Do not write on, or poke the product with sharp pen tips.

Before cleaning, make sure that the product has been turned off and that the adapter has been pulled out from the outlet. Clean with a dry cloth.

Keep the product clean and store in a safe place. Keep away from environments with high humidity or temperature.

MALFUNCTIONING AND MAINTENANCE

- If the product cannot be started:
 - (1) Fully charge the product, which takes 4 hours.
- If the product is unable to charge:
 - (1) Reconnect the charger to the power port;
 - (2) Try a different power socket.
- If the product stops during use:
 - (1) Auto timer has gone off;
 - (2) It has stopped due to excess pressure. Let it rest before restarting.
- If the product appears damaged:
 - (1) Do not continue to use. Consult dealer or technical personnel with maintenance experience for help.

NOTE: Do not disassemble or modify the product by yourself. NAIPO will not be responsible for any loss caused by unauthorised disassembly or modification of the product. This action may put the user in danger and/or render the product unusable. Please contact NAIPO customer service for further assistance.

NAIPO shall not be liable for costs, damages or repairs incurred as a result of, due to, or from:

- The warranty period having expired.
- Unauthorized actions such as: product repair, misuse, abuse, modification, disassembly, etc.
- Artificial or accidental damage.
- Damaged internal and external packaging or missing accessories.
- Damage caused by force majeure.

This product has an internal, non-replaceable battery. This battery is not user-serviceable. Please dispose in accordance with local, state, province, and country regulations.

FCC COMPLIANCE

This device complies with Part 15 of the FCC Rules. Operation is subject to the following two conditions:

- (1) This device may not cause harmful interference, and
- (2) This device must accept any interference received, including interference that may cause undesired operation.

The manufacturer is not responsible for any radio or TV interference caused by unauthorized modifications to this equipment. Such modifications could void the user authority to operate the equipment.

NOTE: This equipment has been tested and found to comply with the limits for Class B digital device, pursuant to part 15 of the FCC Rules. These limits are designed to provide reasonable protection against harmful interference in a residential installation. This equipment generates, uses and can radiate radio frequency energy and, if not installed and used in accordance with the instructions, may cause harmful interference to radio or television reception, which can be determined by turning the equipment off and on, the user is encouraged to try to correct the interference by one or more of the following measures:

- Reorient or relocate the receiving antenna.
- Increase the separation between the equipment and the receiver.
- Connect the equipment into an outlet on a circuit different from that to which the receiver is connected.
- Consult the dealer or an experienced radio/TV technician for help.

Compliance information statement shall issue by a US Local Responsible Party. Caution: Changes or modifications not expressly approved by the party responsible for compliance could void the user's authority to operate the equipment.

SPECIFICATIONS

PRODUCT NAME	oPillow Plus Back Massager
MODEL	oPillow-P1
NET WEIGHT	1.7kg 3.7lb
DIMENSION	L335*W200*H110mm 13.0*8.0*4.0in
BATTERY	2600mAh
INPUT	12V==2.0A
RATED POWER	24W

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